

CAPP Wash™



Easy as never before

Due to CAPP's special two-in-one design of the washing tubes, the CappWash makes hand-held washing very easy: Just lower the tubes into the wells, press the knob for the desired number of washings and move on to the next row. It's as simple as that!

Have your own washer

Since the CappWash is very compact and budget friendly, you now have the possibility of having a washer entirely to yourself. No waiting time because the department washer is occupied or down, no programming to be done before you can start. The simple and efficient design makes the CappWash almost completely free of mechanical problems - and with no electrical connections, there will never be software problems...!

Ideal for your small scale work

Whether you work in research or clinical laboratories, we believe that the CappWash is the ideal washing tool for your small scale immuno assays. Even if you work with large scale immunological testing you may wish to have a CappWash in reserve, just for the case that the automated washer is out of work - or for your pilot series and test applications.

Made for you!

- Easy to use
- Very robust and reliable
- Uniform and efficient washing
- Less wash solution residuals
- Wash strips, 96-well and 384-well plates
- Adjustable to all well shapes
- Autoclavable
- Budget friendly



Two-in-one washing tubes



Screws for easy well-depth adjustment




Push the button to fill - release to empty



A holder follows with each CappWash

Cat. no.	Description
W-8	8 Channel Washer for 96-well microplates or strips
W-12	12 Channel Washer for 96-well microplates or strips
W-16	16 Channel Washer for 384-well microplates
WP-1	Vacuum pump for CappWash 8, 12, and 16
WB-1	Bottle for washing solution
WB-2	Vacuum trap / waste collection bottle

Connection chart

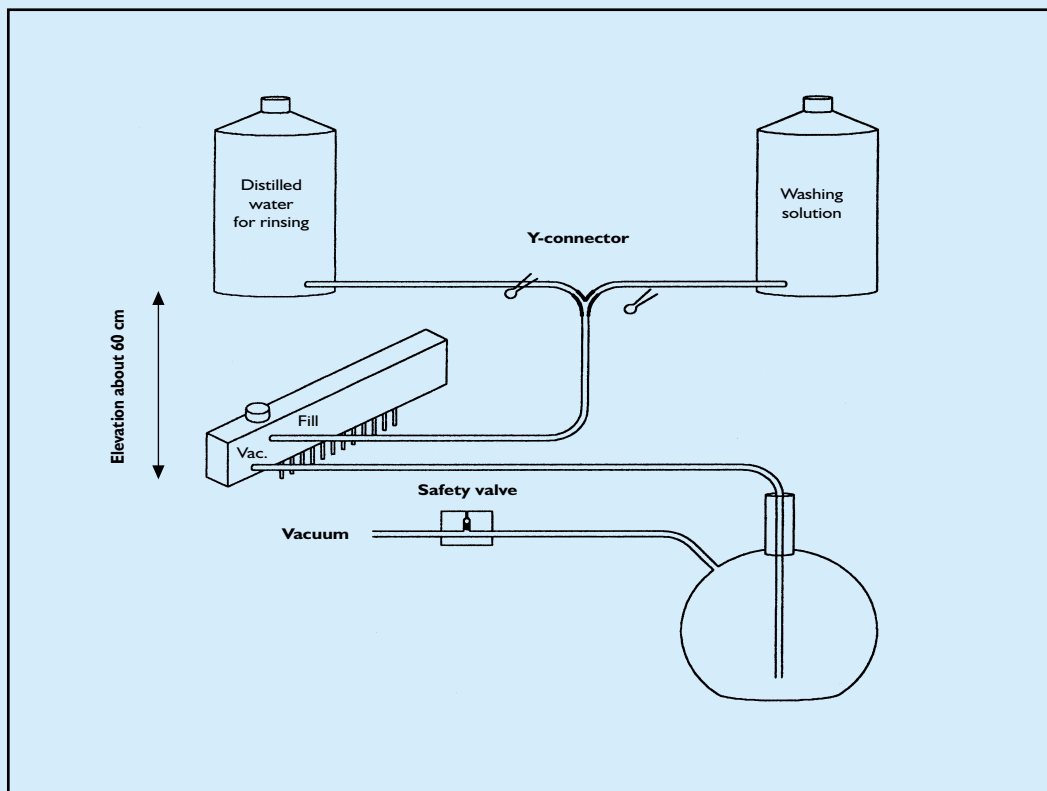

Distributor



Cappelen Laboratory Technics A/S